

# Retrospective PeriCALM Patterns™

## Installation Guide


- Double-click on the **Retrospective PeriCALM Patterns Setup** icon to install the application.



- Once the installation is completed, click **Finish** and Retrospective PeriCALM Patterns will launch automatically. Alternately, you can launch the application from the desktop icon shown below, or by clicking on **Start → All Programs → Retrospective PeriCALM Patterns**.



- Before using Retrospective PeriCALM Patterns for the first time, you must activate the application. Click on **Request Activation Code**.

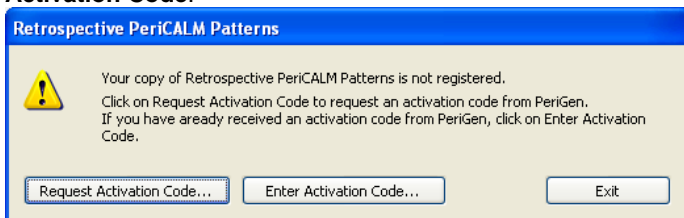

- Fill in the information requested and then click **Send**. This sends your information to PeriGen via email. Alternately, you can manually copy and paste the information by clicking **Copy to Clipboard** and pasting the contents into an email that should be sent to [patterns@perigen.com](mailto:patterns@perigen.com).

A dialog box titled 'Request Activation Code' with a blue header. The text inside says: 'Send the following information to patterns@perigen.com'. There are four input fields: 'Name:', 'Organization:', 'Telephone:', and 'Serial ID:'. The 'Serial ID' field contains the text '3A-AE-AF-54-AC-FB-EC-DA'. Below the fields are two buttons: 'Send' and 'Copy to Clipboard', separated by the word 'or'. At the bottom, there is a 'Close' button. A note at the bottom says: 'Once you have sent the email, you may close this window. When you receive your activation code, run Retrospective PeriCALM Patterns again and select "Enter Activation Code".'

- Within a short period, PeriGen will email you an activation code which is specific for your computer. Open Retrospective PeriCALM Patterns once again, this time clicking on **Enter Activation Code**. Enter your activation code and click **OK**.

A dialog box titled 'Enter Activation Code' with a blue header. It has a single input field labeled 'Activation Code:'. At the bottom are two buttons: 'OK' and 'Cancel'.

- Retrospective PeriCALM Patterns is now activated and will remain activated for the period of time agreed on between you and your PeriGen sales rep.
- To learn more about *Retrospective PeriCALM Patterns* or to learn more about the real-time version of this application **PeriCALM Patterns 1.14**, click on the icon, or contact us at [sales@perigen.com](mailto:sales@perigen.com) or 1-877-263-2256.