


Retrospective PeriCALM Patterns™

Quick Start Guide

Opening a Tracing File


- Click **Open File** on the Welcome Screen

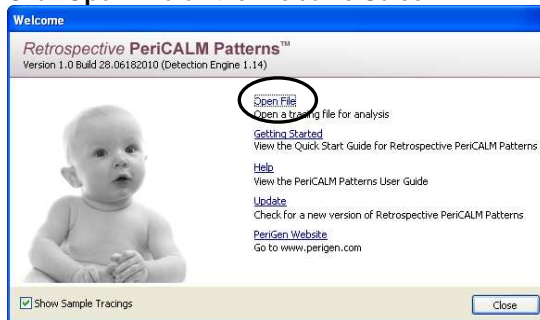

- If the Welcome Screen is already closed or if you wish to open additional tracing files, click **Open** in the bottom toolbar.




- Select your **.v01** file and click **Open**. Please wait as your file is analyzed for baselines, contractions, accelerations and decelerations. *Note: For very long tracings, only the most recent 48 hours are displayed.*

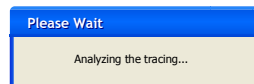

Understanding the Display

- The **Patient Banner** ⁽¹⁾ displays the filename of the tracing currently being displayed. *Note: Retrospective PeriCALM Patterns neither reads nor stores any clinical or demographic patient data.*
- The **Patient List** ⁽²⁾ shows all tracing files currently open.
- The tracing is displayed as both an **Expanded View** ⁽³⁾ showing 17 minutes of tracing and a **Compressed View** ⁽⁴⁾ showing between 2 and 12 hours of tracing.
- Contractions** are represented as beige lines, **Accelerations** as green lines and **Decelerations** as yellow lines. Clicking on any of these lines yields additional information.
- The **Persistent Contractility Index** ⁽⁵⁾ displays the patient's level of contractility and indicates when it is excessive (as defined by your facility).
- The **Slider Window** ⁽⁶⁾ highlights a specific 17 minutes of tracing and shows mean contraction interval, mean baseline and mean baseline variability for that time period.
- All *Retrospective PeriCALM Patterns* features are accessed via the **Toolbar** ⁽⁷⁾.

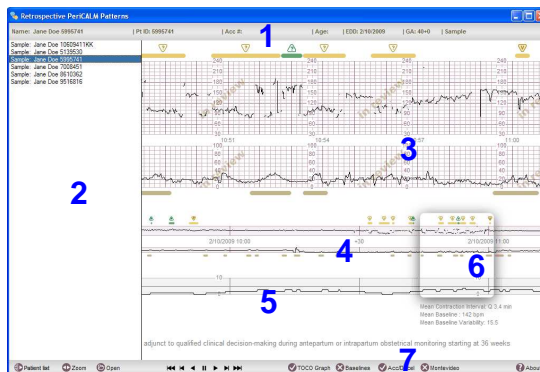

Learn More

- To learn more about *Retrospective PeriCALM Patterns* or to learn more about the real-time version of this application **PeriCALM Patterns 1.14**, click on the **About** icon, or contact us at sales@perigen.com or 1-877-263-2256.